

BATRIK

ARC™ WIRE STYLE BRUSHES

SINGLE-USE AND REUSABLE

You've got options with our ARC™ Wire Brushes. Our single-use and reusable brushes are available in both flexible and rigid varieties and come in a multitude of sizes. ARC™ Wire Brushes are guaranteed to last longer than twisted wire and preserve the integrity of your valuable instruments.

BRUSH SELECTION TIPS:

INSTRUMENT SPECIFICATIONS

What is the maximum length of the instrument(s) to be cleaned?

Add at least 2" (5cm) to that length, so you have an adequate length to work with.

What is the inner diameter of the instrument?

Add at most 0.125" (3mm) to determine the diameter of the bristle to ensure effective cleaning. The bristle diameter has to slightly exceed the inner tube diameter. Bristles that

are too large will bend too much and not clean effectively.

BRUSH SELECTION

Simply find the section in our catalogue that best corresponds to the overall length and diameter of your instrument. As a rule, if your present brush is 8" (20.32cm) in length and we only offer 10" (25.4cm) length with your required diameter, that is entirely acceptable.

In many cases, our ARC™ Wire Brushes will have more bristles on the brush. As an example, many 5mm (0.19") diameter brushes are only 1" (2.5cm) in length. At Batrik, we strived to provide the greatest amount of bristles to ensure the most effective cleaning.

CUSTOMIZABLE WIRE BRUSHES


A unique aspect of single-use ARC™ Wire Brushes is that they can be customized to any length:

- Up to 2000mm for Wire Brushes on plastic semi-rigid sheaths, rigid brushes from 5 to 20mm in diameter.
- No limit in length for Wire Brushes on flexible plastic sheaths, flexible brushes from 2-5mm in diameter

RIGID



FLEXIBLE

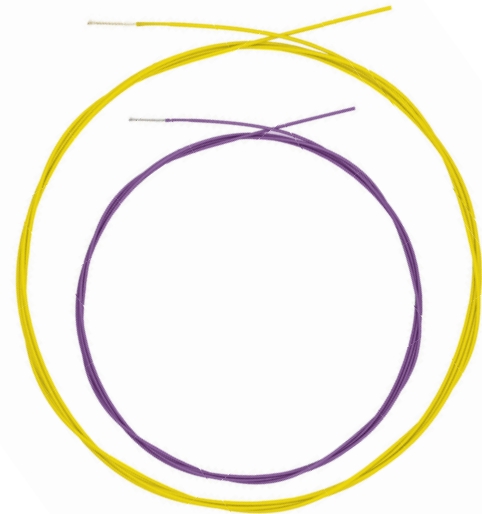

12" / 30 CM SINGLE-END WIRE BRUSH

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-31502-S	41-31502	White	12.401" / 31.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	41-32503-S	41-32503	Yellow	12.795" / 32.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	411-32504-S	411-32504	Violet	12.795" / 32.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	411-32505-S	411-32505	Grey	12.795" / 32.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
RIGID BRUSH & SHEATH	41-35005-S	41-35005	White	13.779" / 35 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
	41-35006-S	41-35006	Green	13.779" / 35 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	41-35008-S	41-35008	Blue	13.779" / 35 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-35010-S	41-35010	Red	13.779" / 35 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-35012-S	41-35012	Yellow	13.779" / 35 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-35015-S	41-35015	Green	13.779" / 35 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-35020-S	41-35020	Blue	13.779" / 35 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING Single-Use: Bag of 20 / Reusable: Bag of 3

18" / 45 CM SINGLE-END WIRE BRUSH

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-46502-S	41-46502	White	18.307" / 46.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	41-47503-S	41-47503	Yellow	18.7" / 47.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	411-47504-S	411-47504	Violet	18.7" / 47.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	411-47505-S	411-47505	Grey	18.7" / 47.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
RIGID BRUSH & SHEATH	41-50005-S	41-50005	White	19.685" / 50 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
	41-50006-S	41-50006	Green	19.685" / 50 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	41-50008-S	41-50008	Blue	19.685" / 50 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-50010-S	41-50010	Red	19.685" / 50 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-50012-S	41-50012	Yellow	19.685" / 50 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-50015-S	41-50015	Green	19.685" / 50 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-50020-S	41-50020	Blue	19.685" / 50 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Single -Use: Bag of 20 / Reusable: Bag of 3

24" / 60 CM SINGLE-END WIRE BRUSH



	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-61502-S	41-61502	White	24.212" / 61.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1mm
	41-62503-S	41-62503	Yellow	24.6" / 62.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1mm
	411-62504-S	411-62504	Violet	24.6" / 62.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55mm
	411-62505-S	411-62505	Grey	24.6" / 62.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75mm
RIGID BRUSH & SHEATH	41-65005-S	41-65005	White	25.59" / 65 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3mm
	41-65006-S	41-65006	Green	25.59" / 65 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3mm
	41-65008-S	41-65008	Blue	25.59" / 65 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-65010-S	41-65010	Red	25.59" / 65 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-65012-S	41-65012	Yellow	25.59" / 65 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-65015-S	41-65015	Green	25.59" / 65 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-65020-S	41-65020	Blue	25.59" / 65 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Single-Use: Bag of 20 / Reusable: Bag of 3

30" / 75 CM SINGLE-END WIRE BRUSH


	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-76502-S	41-76502	White	30.118" / 76.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	41-77503-S	41-77503	Yellow	30.511" / 77.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	411-77504-S	411-77504	Violet	30.511" / 77.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	411-77505-S	411-77505	Grey	30.511" / 77.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
RIGID BRUSH & SHEATH	41-80005-S	41-80005	White	31.496" / 80 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
	41-80006-S	41-80006	Green	31.496" / 80 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	41-80008-S	41-80008	Blue	31.496" / 80 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-80010-S	41-80010	Red	31.496" / 80 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-80012-S	41-80012	Yellow	31.496" / 80 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-80015-S	41-80015	Green	31.496" / 80 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-80020-S	41-80020	Blue	31.496" / 80 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Single-Use: Bag of 20 / Reusable: Bag of 3

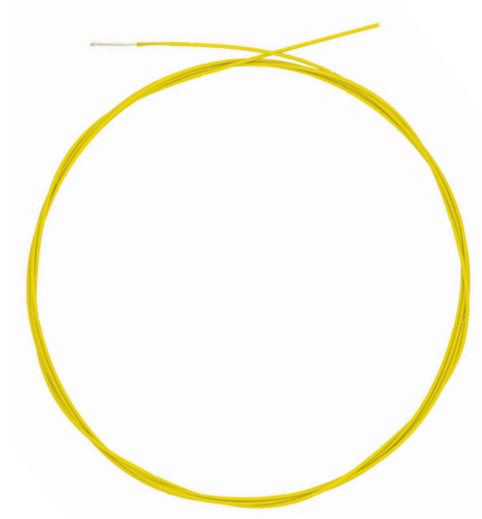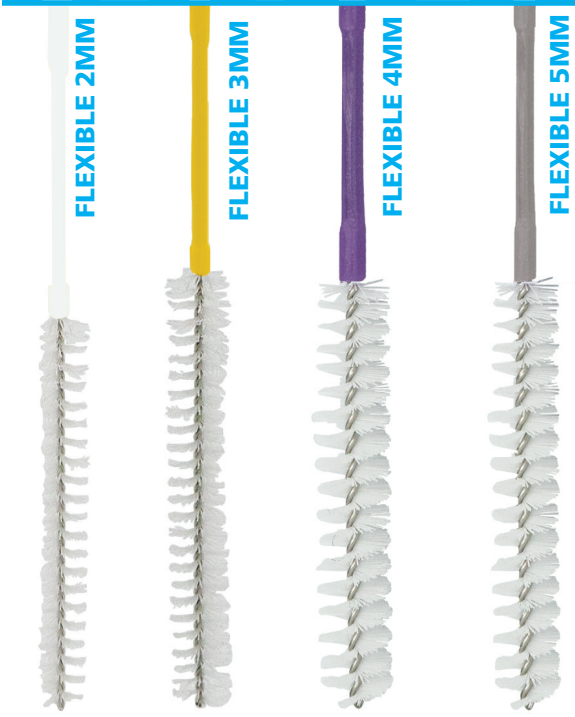

FLEXIBLE

Flexible brush diameter (Ø) ranging from 2 mm to 5 mm

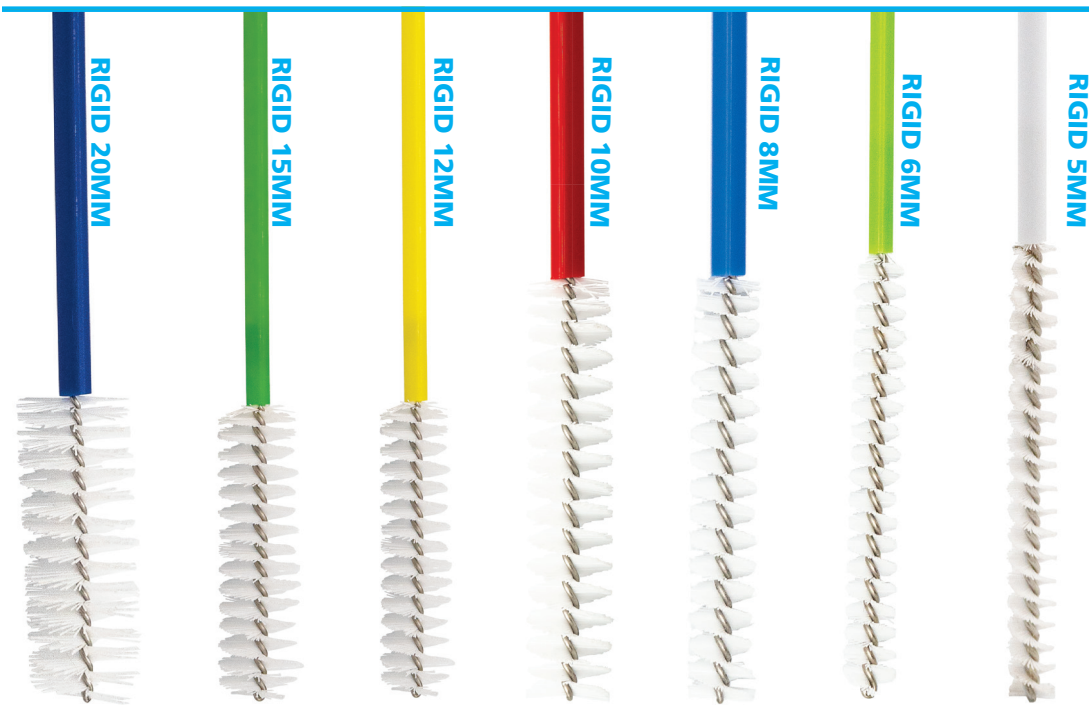

RIGID

Rigid brush diameter (Ø) ranging from 5 mm to 20 mm

BRUSH SELECTION SYSTEM (Actual Size)

2mm 0.078	3 mm 0.118"	4 mm 0.157"	5 mm 0.196"	6 mm 0.236"	8 mm 0.314"	10 mm 0.394"	12 mm 0.472"	15 mm 0.591"	20 mm 0.787"